

Frozen in Time Ltd.

Manufacturers of Freeze Drying Machines
and Vacuum Cold traps

Lablyo Extra Freezedrier

UK Manufactured

The **LablyoExtra** offers capabilities similar to our larger H series pilot freeze driers in a small bench top system.

These include:

Control panel for programming steps using time and temperature.


Data available to download onto a memory stick.

3 temperature controlled stainless steel shelves with circulation fluid eliminates the need to pre freeze samples and allows temperature control during the freeze drying process.

4 mobile temperature probes in drying chamber for product sampling.

Condenser capacity 3kg, temperatures to -80°C

Pneumatic shelf stoppering for vial closing.

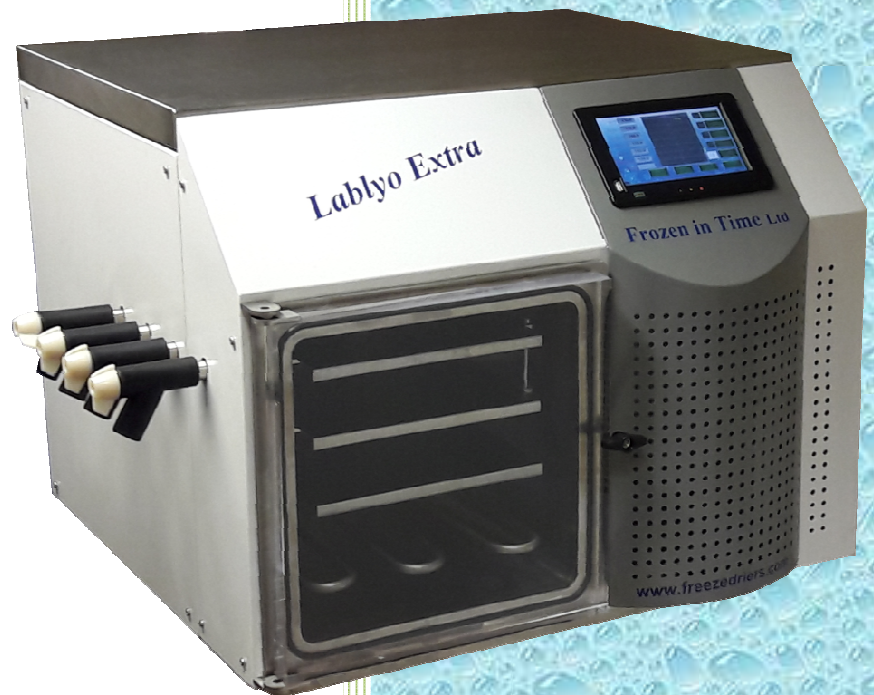

External 4 port manifold to allow for flask freeze drying.

Application:

Laboratory research and development freeze drying.

96 well plates.

Serum bottles or vials.

Bulk products for analysis.

Small R&D freeze drying.

The **Lablyo Extra** comes with clear acrylic door for observing the process.

Overall Dimensions: Height 55cm x length 70cm x width x 85cm

The **Lablyo Extra** includes a vacuum pump and connection hose that sits separately.

FROZEN IN TIME Ltd

York Road Industrial Park, Sheriff Hutton, York, YO60 6RZ England.
TEL.01347 878158.E-MAIL.info@freezedriers.com WEBSITE. www.freezedriers.com